





stretcher series  
**P-172**  
 basic transport stretcher

The P-172 Series general transport stretcher provides an economical patient transport solution. This non-hydraulic stretcher series features a 1<sup>3</sup>/<sub>8</sub>" heavy-gauge O.D. tubular frame, continuous rubber bumper, long-life, double ball bearing, 8" swivel casters for quiet movement, and convenient hand holds on each corner.

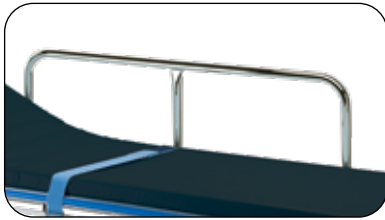


**RAIL SYSTEM**

Easy to wipe clean and disinfect side rails.

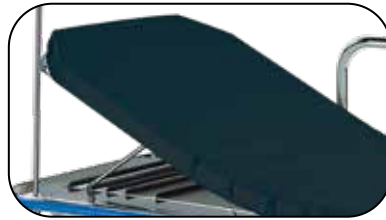

**FOWLER HEAD REST**


Fowler provides additional assistance for the safe and easy raising and lowering of the patient.

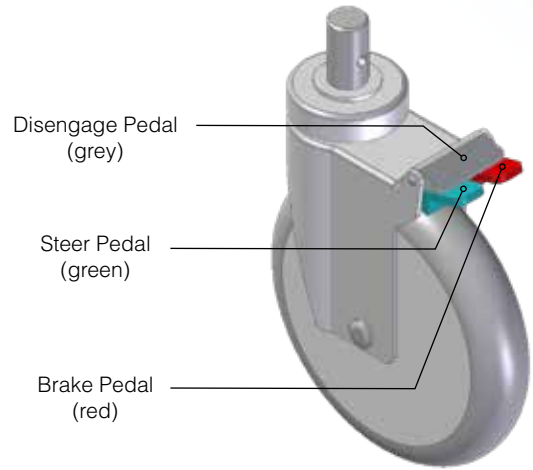

**CASTERS**



I.V. Poles Socket Locations

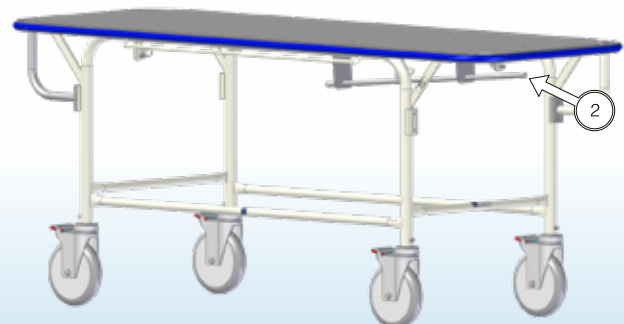

I.V. Pole Storage Tube Location

## STANDARD FEATURES

- 350 lbs. weight capacity
- Designed as a cost efficient transport stretcher
- Four Brake & Steer Casters
- Heavy-gauge, 1<sup>3</sup>/<sub>8</sub>" O.D. tubular frame
- Fixed litter top height
- Full perimeter bumper
- Convenient hand holds on each corner
- Long-life polyurethane, double ball bearing, 8" swivel casters for quiet movement
- Latex-free



## SPECIFICATIONS

Model Number	P-172
Weight Capacity	350 lbs. (159 kg)
Overall Length	78" (198 cm)
Overall Width	28" (71 cm)
Fixed Height	33 <sup>1</sup> / <sub>4</sub> " (85 cm)
<b>Optional Litter Positioning</b>	
Fowler Backrest	3°, 12°, 26°, 38°, 44°
Patient Surface	77" x 26" (196 x 66 cm)
Optional Side Rails	12" x 43" (30 x 109 cm)
Casters	8" (20 cm)
Warranty	1 year parts and labor

## P-172 CONFIGURATIONS

	Swivel Caster	Brake & Steer Caster	P-170-1 Shelf	P-170-2 Rails	P-170-3 Fowler	P-170-4 IV Pole	P-170-5 Strap	P-170-6 Bumpers	P-170-8 Mattress	P-170-9 O <sub>2</sub> Tank Holder
P-172		4								
P-172-C		4	•	•	•	•	•		•	•

## OPTIONS

### Accessory

#### Storage Shelf:

Shelf ..... P-170-1

#### Side Rail:

Stainless Steel..... P-170-2

#### Fowler:

Stainless Steel..... P-170-3

### Accessory

#### IV Pole:

Fixed Height ..... P-170-4

#### Restraint Strap:

Velcro..... P-170-5

#### Rolling Bumpers:

7<sup>1</sup>/<sub>2</sub>" round rubber..... P-170-6

### Accessory

#### Mattress:

Pad (2<sup>1</sup>/<sub>2</sub>" thick black ) . P-170-8

#### O<sub>2</sub> Holder:

E Size Tank Holder ..... P-170-9

(requires storage shelf P-170-1)

Bed Exit Monitor ..... 750216




### MATTRESS PAD

Soft-vinyl 2<sup>1</sup>/<sub>2</sub>" thick black mattress pad measures 26"W x 77"L.



### O<sub>2</sub> HOLDER

Bracket for E-size oxygen tank. Requires lower shelf.

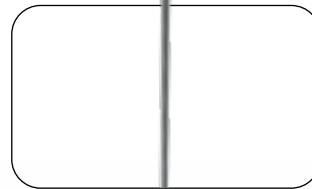


### IV POLE

Two hook pole includes holder and storage tube.



### FWOWER HEAD REST

Stainless steel construction with an adjustable range from 0° - 45°.

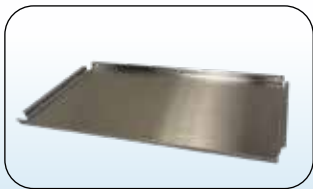

### LOWER SHELF

Stainless steel storage shelf.




### RESTRAINT STRAP

Blue Velcro restraint strap.



### CASTERS

8" double ball bearing casters.

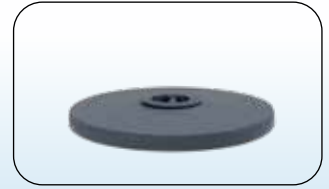

### BUMPERS

7<sup>1</sup>/<sub>2</sub>" rolling rubber bumpers