




stretcher series  
**5400**  
 universal procedure stretcher

The 5400 Series universal transport stretcher is versatile, yet economical. It's the right stretcher for ED, OR, PACU, or any department where cost, efficiency and quality are required.



**RAIL SYSTEM**

Easy to wipe clean and disinfect, full length tuck-away side rails allow for easier patient transfer by creating a zero transfer gap. Rails feature a release latch centrally located on either side of the litter top under the patient surface.

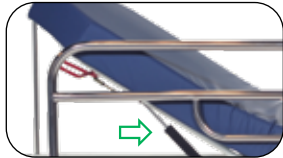

**FOWLER ASSIST**

Pneumatic assist fowler provides additional assistance for the safe and easy raising and lowering of the patient.

**Pneumatic Assist Fowler**

400 N (90#) force



**ZERO TRANSFER GAP**

Tuck-Away side rails allows the caregiver to position the 5400 Series Stretchers flush with the bed for a zero transfer gap and easier patient transfer. No expensive or cumbersome transfer devices are required.



**HYDRAULIC PEDALS**

Hydraulic Foot Pedals are located on each side of the stretcher within the patient care zone. Color-coded pedals and function graphics are quickly and easily understood.




**COLUMN SUPPORT**

Steel center column support with single hydraulic lift cylinder design uses sixteen precision ball-bearings to smoothly raise and lower the patient.

**PATIENT LIFTS**

The 5400 Series Stretcher has 5<sup>7</sup>/<sub>8</sub>"H" of clearance between floor and cowling to easily accommodate and utilize patient lifts.

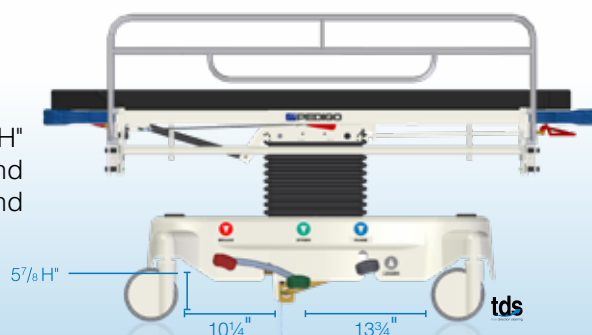

**TRUE DIRECTION STEERING**

For improved stretcher control engage the True Direction Steering system. The steer caster located at the patient's right foot end of the stretcher will lock once the steering pedal is applied.




## STANDARD FEATURES

- 500 lbs. weight capacity
- True Direction Steering (TDS)
- Side rails - full length, single release
- Center column design with precision ball bearings
- Single Trendelenburg assist cylinder (positive/negative - pneumatic)
- Pneumatic assist fowler 400 N (90#) force
- Dual side-mounted hydraulic foot pedals
- Brake & Steer activated foot pedals
- Four corner IV pole receptacles
- Four restraint strap positions
- High impact corner bumpers
- 8" easy mobility casters
- Latex-free

## SPECIFICATIONS


Model Number	5400 <sup>N</sup>	5400 <sup>W</sup>
Weight Capacity	500 lbs. (227 kg)	500 lbs. (227 kg)
Overall Length	83" (211 cm)	83" (211 cm)
Overall Width	29.5" (75 cm)	33.5" (85 cm)
<b>Height Range</b>		
Low	23" (58 cm)	23" (58 cm)
High	32.75" (83 cm)	32.75" (83 cm)
<b>Litter Positioning</b>		
Fowler Backrest	0°-90°	0°-90°
Trendelenburg	+12°	+12°
Reverse Trendelenburg	-7°	-7°
Patient Surface	77" x 27" (196 x 69 cm)	77" x 31" (196 x 77 cm)
Side Rails	15" x 63 <sup>5</sup> / <sub>8</sub> " (38 x 162 cm)	15" x 63 <sup>5</sup> / <sub>8</sub> " (38 x 162 cm)
Casters	8" (20 cm)	8" (20 cm)
Warranty	1 year parts and labor	1 year parts and labor

## OPTIONS

Accessory	N	W	Accessory	N	W	Accessory	N	W
<b>X-Ray Board</b>			<b>Side Rail Options:</b>			<b>Knee Flex</b>		
Full.....	750203	750204	Rail Pad (pair) .....	5716002	same	Hydraulic.....	750208	750209
Partial .....	750205	750206	<b>IV Pole Options:</b>					
<b>Mattresses:</b>			Fixed Height.....	2101	same	Leg Lift	Knee Flex	Trauma
3" Standard (black).....	5828002	5831002	Adjustable / Removable .....	5984001	same	<b>Additional Options:</b>		
3" Premium (blue).....	5850002	5852002	Adjustable / Attached.....	750501	750502	Bed Exit Monitor .....	750216	same
4" Standard (black).....	5829002	5832002	Fold Down .....	59175001	same	Extra Release Pedal - Hydraulic for stretchers with <b>TDS</b> ....	750105	same
4" Premium (blue).....	5851002	5853002	IV Pole Caddy .....	59163001	same	Storage Bag Hook.....	750215	same
4" Pressure Relief .....	5833002	5834002	<b>Push/Pull Bar Options:</b>			Paper Roll Holder .....	5966001	5966003
5" Premium (blue)† .....	5855002	5856002	Push/Pull Handles (CuNi alloy).....	750223	750224	Quick Release O <sub>2</sub> Holder .....	750213	same
† <i>Special order only.</i>			Push/Pull Bar (Removable) ..	5985001	5986001	O <sub>2</sub> Tank Holder (Vertical) ..	5981001	same
<b>Restraint Straps Options:</b>			Bar Pad .....	5704001	5705001	Head/Foot Extension .....	5996001	same
Velcro (pair).....	5903001	same	Chart Holder.....	5967001	same	Head/Foot 3" Pad .....	5708001	same
Buckle (pair).....	5904001	same	Monitor Tray Board.....	59115001	59116001	Stirrups Set.....	5910001	same
			<i>(includes chart holder 5967001)</i>					



SIDE RAIL PADS

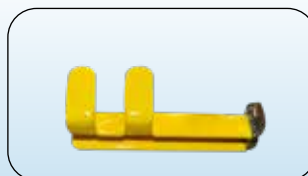

STORAGE BAG HOOK



RESTRAINT STRAPS